



zoomax

by  control-x
MEDICAL

cxmed.com

ZOOMAX GOLD VETELEV SYSTEM



VETELEV ELEVATING TABLE

A SIMPLIFIED AND COST-CONSCIOUS VERSION OF OUR POPULAR PHOENIX MEDICAL PATIENT TABLES, FEATURING



SUBSTANTIAL FREEDOM OF POSITIONING

- 36 cm vertical movement
- Floating table top with electromagnetic brakes (152 cm x 71 cm with extension option)
- Over 100 kg lifting capacity




TOTAL SAFETY AND COMFORT

- Integrated crash guard sensors on the table
- Grid cabinet for standard cassette size wifi or fixed DR panels ensuring crystal clear images
- Configuration for fixed or removable grids, including safe „parking slot” for unused grid
- Wireless foot pedal (with long-life battery) to minimize cabling (optional)
- Foot pedal doubles as light field switch



ZOOMAX TUBE STAND

- Sleek, robust design
- Fixed position tube with 100-136 cm SID (with Vetelev table)
- Free-standing setup for total positioning freedom (optional)

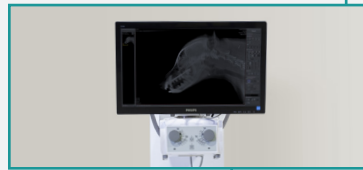

SELECTION OF CONTROL CONSOLES

- 22" Touchscreen Monitor Console which can be used to set exposure parameters and for full image capture and analysis
- ZooMax LC or the CPI Rad Console which are membrane consoles offering the comfort of all key generator parameter settings, including APRs and advanced controls such as AEC




INTEGRATED GENERATOR CONTROL AND IMAGE ACQUISITION

Single console seamless workflow:


- Unlimited APRs
- Full DICOM compatibility

Sophisticated image processing features:

- Detector status indicators
- Wifi signal strength indicator
- Wide variety of DR detector choices
- Portable solutions with tablet, laptop workstation
- Mini PACS included
- Cloud storage available
- Mobile viewing on any device

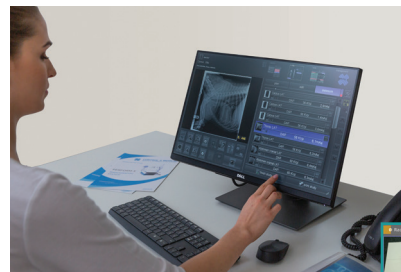

Patient positioning guide

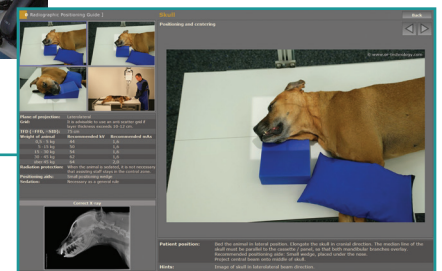

HIGH QUALITY COMPONENTS



- Industry-leading x-ray tubes
- High-density grids
- Certified manual collimators with light beam or laser for easy positioning
- Powerful generators for minimal exposure times

GENERATOR OPTIONS

- 32 kW to 50 kW output range
- Stored energy generators up to 50 kW (capacitor discharge) at 15 Amp
- 2.8 kW to 5 kW monoblock generators

REMOTE SUPPORT

Service software allows troubleshooting to minimize downtime.

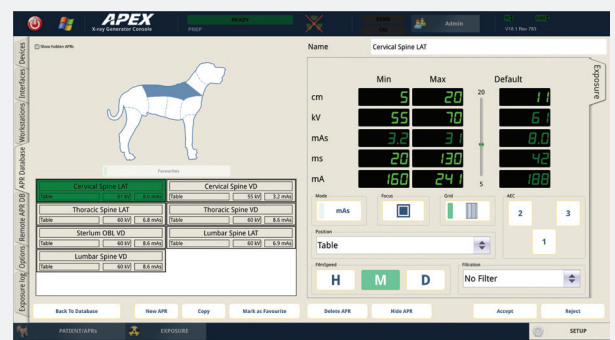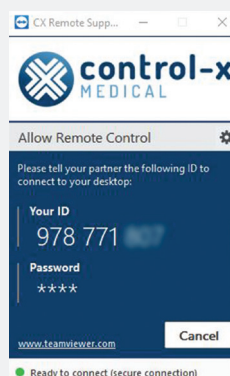

Please contact us at sales@cxmed.hu or sales@cxmed.com

CERTIFICATIONS

Control-X Medical veterinary diagnostic products are CE marked.

